


Askina® Foam Range

Hydrophilic Foam Dressings

B. Braun Wound Care. Right. From the start.

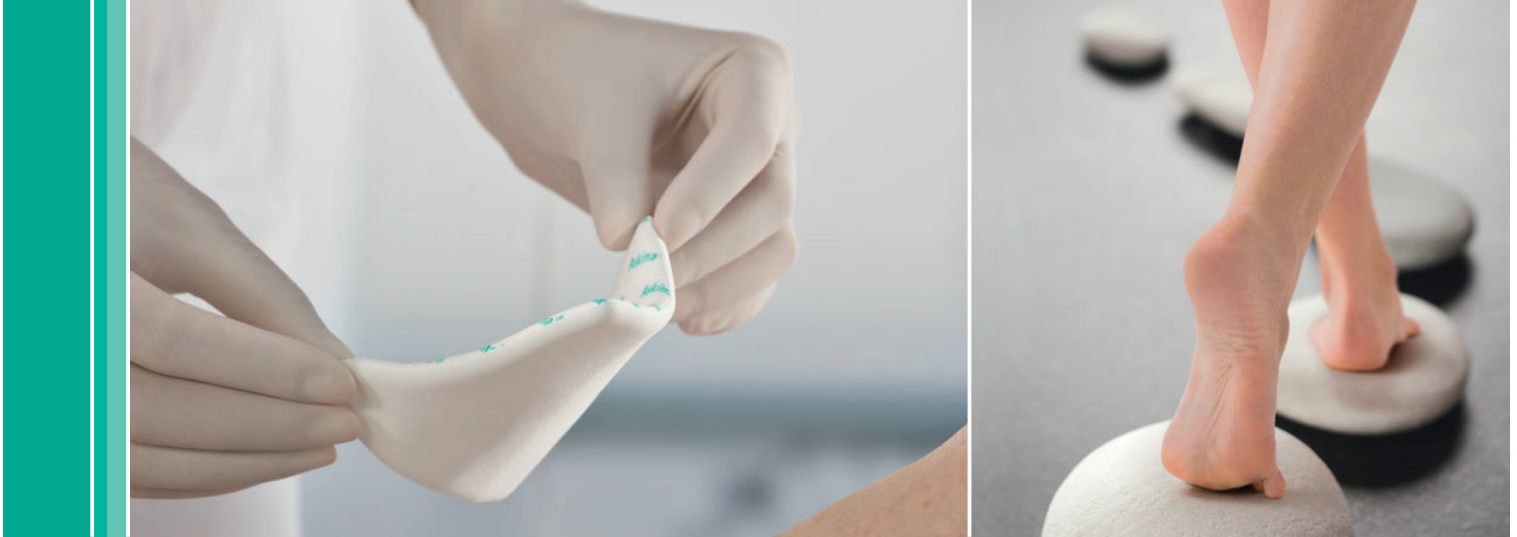

Wound Management

B | BRAUN
SHARING EXPERTISE



Established 1922

Medical and Surgical Requisites Pty Ltd

Call within Brisbane 07 3859 2900
Phone 1300medsurg (1300 633 787)

33 Fulcrum Street Richlands Queensland 4077
Email medical@medsurg.com.au

www.medsurg.com.au
Fax 07 3859 2995



Askina® Foam | Askina® Foam Cavity | Askina® Cavity Strips

Hydrophilic foam dressings



Excellent fluid handling characteristics

Cost efficient: no need for frequent dressing changes

Enhances a moist wound environment

No adherence to the wound;
atraumatic dressing changes

Visual control of exudate absorption

Does not leave pressure marks when used under
compression therapy

Askina® Foam Cavity is highly resistant:
easy to remove, without leaving particles

Askina® Foam


is a two layered non-adherent foam dressing consisting of

- › A soft hydrophilic polyurethane foam layer which is breathable and has a high fluid handling capacity (199 g / 100 cm²)
- › A thin, transparent and protective polyurethane film, which is water impermeable and bacteria resistant

Askina® Foam Cavity and Askina® Cavity Strips

are made of the same foam material as Askina® Foam, without the polyurethane backing. Their specific shape, high absorption capacity and resistance make them ideal for dressing deep, highly exuding wounds as well as small wound sinuses (Askina® Cavity Strips).

Innovative product design

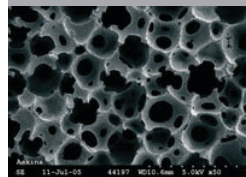

The specific open cell structure of Askina® Foam allows rapid and vertical absorption of exudate, with no risk of in-growth of newly formed granulation tissue.



Established 1922

Medical and Surgical Requisites Pty Ltd

Call within Brisbane 07 3859 2900 33 Fulcrum Street Richlands Queensland 4077 www.medsurg.com.au
Phone 1300medsurg (1300 633 787) Email medical@medsurg.com.au Fax 07 3859 2995



Indications Askina® Foam

Moderately to heavily exuding wounds

- › Stage I–IV pressure ulcers
- › Venous and arterial leg ulcers
- › Diabetic foot ulcers
- › 1st and 2nd degree burns
- › Traumatic wounds

Indications Askina® Foam Cavity


- › Deep cavity wounds with high exudate levels

Indications Askina® Cavity Strips

- › Small cavity wounds, moderately to highly exuding
- › Wound sinuses



- ① Askina® Foam used on a heavily exuding flat wound: visual control of exudate absorption.
- ② Large quantity of exudate is absorbed: vertical absorption, without maceration of the wound edges.
- ③ ④ Combined use of Askina® Foam Cavity, used to fill a deep highly exuding wound, and Askina® Foam, used as an additional absorption dressing.

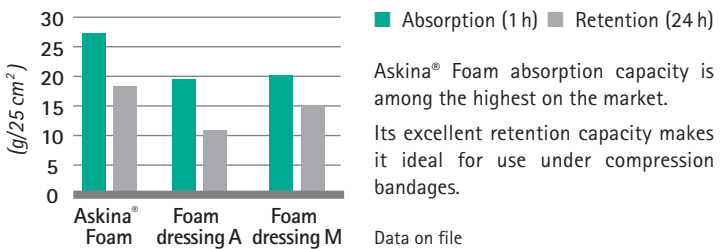

Askina® Foam	Size	Pcs/Pack	Article number
	5 cm x 7 cm	10	7240710
	10 cm x 10 cm	10	7241010
	10 cm x 20 cm	10	7241210
	20 cm x 20 cm	5	7242005
Cavity	2.5 cm x 40 cm	10	7244010
Askina® Cavity Strips	2.5 cm x 40 cm	10	7246010

TIP

- › Askina® Foam can be fixed with tapes or bandages, or with compression bandages
- › Askina® Foam can also be used as a secondary dressing for highly exuding wounds, in combination with Askina® Sorb, or as a secondary dressing for deep wounds, filled with Askina® Sorb Rope, or Askina® Foam Cavity



Medical and Surgical Requisites Pty Ltd

Call within Brisbane 07 3859 2900 33 Fulcrum Street Richlands Queensland 4077 www.medsurg.com.au
 Phone 1300medsurg (1300 633 787) Email medical@medsurg.com.au Fax 07 3859 2995



Askina® Foam Border | Askina® Foam Sacrum | Askina® Trachea

Hydrophilic foam dressings

Askina® Foam Border and **Askina® Foam Sacrum** are made of the same foam material as Askina® Foam with a polyurethane film outer layer, which extends beyond the perimeter of the foam to form an adhesive border.

Askina® Foam Border is a square, two layered foam dressing, particularly adapted for difficult-to-dress areas.

Askina® Foam Sacrum has a triangular shape, well suited for sacral area. It is also suitable for the prevention of skin damages due to incontinence or pressure.

Indications Askina® Foam Border

Moderately exuding wounds

- › Stage I-IV pressure ulcers,
- › Venous and arterial leg ulcers
- › Diabetic foot ulcers
- › 1st and 2nd degree burns
- › Traumatic wounds

Indications Askina® Foam Sacrum

- › Moderately exuding wounds located on the sacral area




Askina® Trachea is a hydrophilic foam dressing with the circular aperture designed to fit neatly around a tracheostomy tube or other drain or stoma.

Indications Askina® Trachea

Moderately to heavily exuding wounds

- › Tracheostomy tubes
- › Drainage sites
- › Stomas

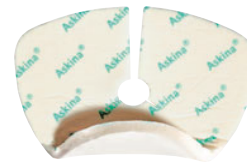

Askina® Foam	Size overall dressing foam island	Pcs/Pack	Article number
Border	16 cm x 16 cm		
	10 cm x 10 cm	10	7231610
	16 cm x 26 cm		
	10 cm x 20 cm	10	7236610
	26 cm x 26 cm		
Sacrum	20 cm x 20 cm	5	7232605
	18 cm x 20 cm		
	14 cm x 16 cm	5	7228205

Askina® Trachea	Size	Pcs/Pack	Article number
	9.5 cm x 8.5 cm	10	7248510



Medical and Surgical Requisites Pty Ltd

Call within Brisbane 07 3859 2900
Phone 1300medsurg (1300 633 787)

33 Fulcrum Street Richlands Queensland 4077
Email medical@medsurg.com.au

www.medsurg.com.au
Fax 07 3859 2995

