

DERM 942

High Frequency Desiccator



The DERM 942 and its patent pending handpiece are designed and engineered to be the most reliable and durable desiccator available today. The unit's life is prolonged by incorporating the latest in chip technology and reducing the number of internal parts, thereby generating less internal heat.

40 watt high frequency desiccator

Adjustable in 1/10th watt increments in the first 10 watts, in monopolar and bipolar

More accurate power with Universal Power Control 100-240 VAC

Less to do and less to buy makes this the most cost effective and user friendly unit on the market

The Bovie DERM 942 is a ground reference unit which means no grounding pad is required.



BOVIE® 942 HIGH FREQUENCY DESICCATOR




STANDARD EQUIPMENT

- A942 (1) electrocautery generator,
- A902 (1) 3-Button Autoclavable Handpiece
- A804 (5) non-sterile sharp dermal tips
- A805 (2) sterile sharp dermal tips
- A806 (5) non-sterile blunt dermal tips
- A807 (2) sterile blunt dermal tips
- A806DE (2) non-sterile DERM-Elite™ Premium blunt dermal tips
- A807DE (2) sterile DERM-Elite Premium blunt dermal tips
- A910ST (2) sterile disposable handpiece sheath
- A910 (2) non-sterile disposable handpiece sheath
- A837 (1) wall mount kit
- 110 VAC Hospital-grade power cord*
- User's guide (MC-55-239-001)*


* Not Pictured



A902

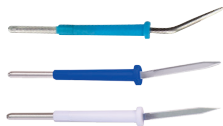

A804, A805, A806, A807, A806DE, A807DE

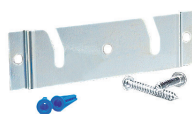

A837



A910ST



A910

COMPLEMENTARY PRODUCTS †



A813
Table Top Stand †



A803
Footswitch †



A827V
Forcep Cable †



Reusable Bipolar
Forceps †



ES20
Ball Electrode †



ES02
Needle Electrode †



A812
Mobile Stand †

† Sold Separately

PRODUCT SPECIFICATIONS

DIMENSIONS

Height: 188 mm
Width: 228 mm
Depth: 105 mm
Weight: < 2.26 kg

OUTPUT POWER

To 40 watts (0-10 watts in
1/10th watt increments on
both monopolar and bipolar)

OUTPUT FREQUENCY

328-492kHz

LINE VOLTAGE

A942: 100-240 VAC

LINE FREQUENCY

50 – 60Hz

WARRANTY

Four Years



Medical and Surgical Requisites Pty Ltd

Call within Brisbane 07 3859 2900 33 Fulcrum Street Richlands Queensland 4077 www.medsurg.com.au
Phone 1300medsurg (1300 633 787) Email medical@medsurg.com.au Fax 07 3859 2995

