


# Bovie

Bovie Medical Corporation's electrosurgery business is now a part of



## BOVIE® BANTAM|PRO™

### *Electrosurgical Generator*

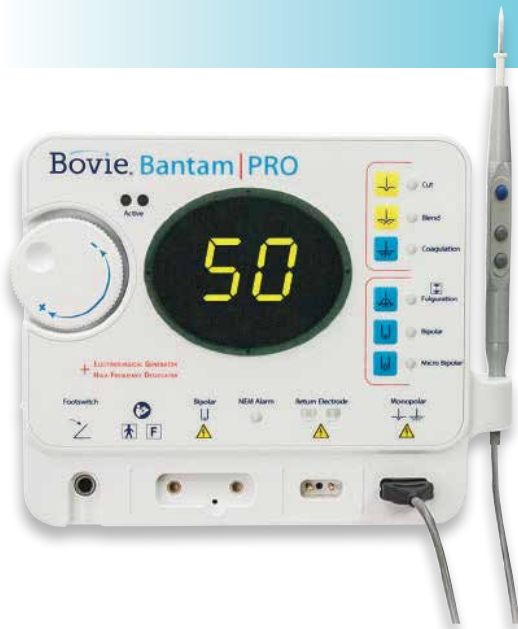

The Bovie Bantam|PRO, A952, is a High Frequency Desiccator with cutting ability. You might call it a split personality. One half of the unit is a high frequency desiccator, an upgrade to our very popular Aaron® 940. The other half is a Bovie Generator, giving you 50 watts of cutting power and coag capabilities – plenty of power for most cutting procedure performed in the doctor's office.

High Frequency Desiccator with 50 watts of Cutting capabilities

Two machines in one convenient package

Adjustable in 1/10th watt increments in the first 10 watts, for the most accurate output available

No grounding pad required when working in the desiccation mode

When in electrosurgical generator mode, the unit features Bovie NEM™ pad sensing technology which monitors the return electrode for optimum patient protection.



**Medical and Surgical Requisites Pty Ltd**

Call within Brisbane 07 3859 2900  
Phone 1300medsurg (1300 633 787)

33 Fulcrum Street Richlands Queensland 4077  
Email [medical@medsurg.com.au](mailto:medical@medsurg.com.au)

[www.medsurg.com.au](http://www.medsurg.com.au)  
Fax 07 3859 2995



# BOVIE® BANTAM|PRO™ A952 ELECTROSURGICAL GENERATOR



## STANDARD EQUIPMENT

- A952 (1) electrosurgical generator,
- A902 (1) handpiece
- A910 (2) non-sterile disposable handpiece sheaths\*
- A910ST (2) sterile disposable handpiece sheaths\*
- ES01 (3), ES02 (1), ES20 (1) sterile disposable electrodes
- A804 (5), A806 (5), A806DE (2) disposable non-sterile dermal tips\*
- A805 (2), A807 (2), A807DE (2) disposable sterile dermal tips\*
- A1252C (2) reusable grounding cord
- ESRE (5) disposable split grounding pads
- A837 wall mount kit\*
- 110 VAC Hospital-grade power cord\*
- User's guide (MC-55-238-001)\*



A902



ES01, ES02, ES20



A1252C



ESRE

\* Not Pictured

## COMPLEMENTARY PRODUCTS †



A813  
Table Top Stand †



A803  
Footswitch †




A827V  
Forcep Cable †



A824  
Reusable Bipolar Forceps †



Lletz Loop  
Electrodes †

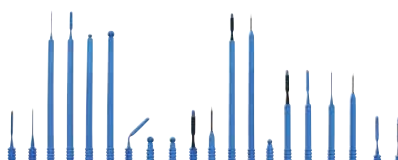

Resistick II  
Coated Electrodes †



Tungsten Wire  
Loop Electrodes †

† Sold Separately

## PRODUCT SPECIFICATIONS

### DIMENSIONS

Height: 188mm (7.4 in.)  
Width: 228mm (8.98 in.)  
Depth: 105mm (4.13 in.)  
Weight: < 2.26 kg (< 5 lbs.)

### OUTPUT POWER

Cut: 50 watts  
Blend: 50 watts  
Coagulation: 50 watts  
Fulguration: 40 watts  
Bipolar: 40 watts  
Micro Bipolar: 40 watts

### OUTPUT FREQUENCY

304 – 451kHz

### LINE VOLTAGE

A952: 100-240 VAC

### LINE FREQUENCY

50Hz

### WARRANTY

Four Years



## Medical and Surgical Requisites Pty Ltd

Call within Brisbane 07 3859 2900 33 Fulcrum Street Richlands Queensland 4077 www.medsurg.com.au  
Phone 1300medsurg (1300 633 787) Email medical@medsurg.com.au Fax 07 3859 2995

