




Easy on-PC

PC spirometer

The **Easy on-PC** is a fully integrated spirometer that leverages our patented TrueFlow technology to deliver precise, accurate, and reliable results.

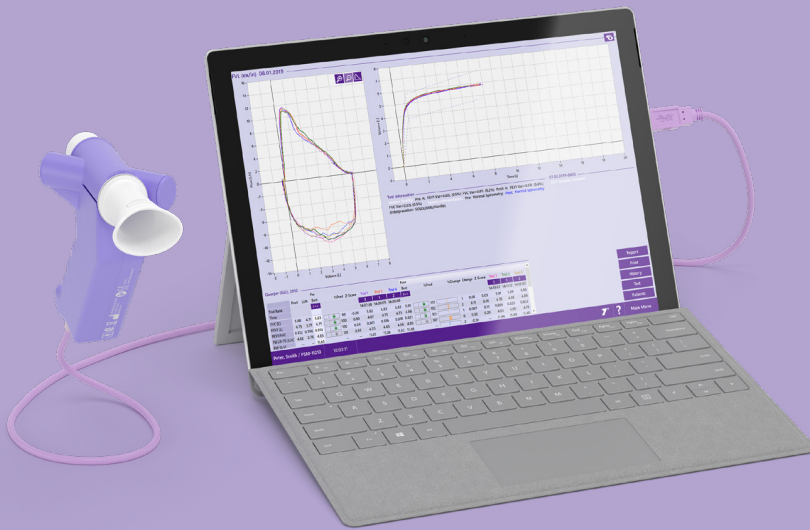

Simple. Certain. Proven.

Features

- Real-time animated incentives
- Proven EMR integration with top EMR/EHR systems
- Long term trending of patient results
- Automated quality feedback and interpretation
- Intuitive software interface
- Integrated provocation test with customizable protocols
- Wide range of selectable predicted values
- Calibration free technology

Incentive screens



Specifications

Test Types: FVC, FVL, SVC, MVV, Challenge and Pre/Post testing

Languages: English, Chinese, Croatian, Danish, Dutch, Finnish, French, German, Italian, Japanese, Norwegian, Portuguese, Russian, Spanish, Swedish, Turkish, Vietnamese

EMR connectivity: HL7, XML, GDT, API with EasyOne Connect

Standards: Meets 2005 & 2019 ATS/ERS standards and NIOSH/OSHA

Hard disk capacity: Installation/system 1 GB. Data up to 4 GB

RAM: 2 GB

EasyOne Connect software compatibility: Windows 7, Windows 8 and 8.1 (32 and 64 bit), Windows 10 (32 and 64 Bit), Windows 11

Part numbers

2700-3: Easy on-PC Spirometer

2700-OCC: Easy on-PC Occupational Health Spirometer Kit
(Sold in the United States only)

2050-1: Spirette box of 50

2050-5: Spirette box of 200

2050-10: Spirette box of 500

TrueFlow
makes the difference

The original ultrasonic flow measurement is highly accurate in all flow ranges, independent of gas composition, pressure, temperature and humidity.

ndd TrueFlow is a hygienic and resistance-free solution that does not require calibration during its lifetime.

EasyOne Connect

EasyOne Connect software offers one database and one platform for all EasyOne point-of-care-solutions. Data management has never been easier.



SCAN ME

ndd Medical Technologies • 300 Brickstone Square, Suite 604 Andover, MA, 01810 USA • Phone: +1-978-470-0923
nnd Medizintechnik AG • Technoparkstrasse 1, CH-8005 Zurich, Switzerland • Phone: +41 44 512 65 00

ID: 9001633 Rev: 01

© Copyright 2022. All rights reserved



Medical and Surgical Requisites Pty Ltd

Call within Brisbane 07 3859 2900 33 Fulcrum Street Richlands Q 4077
Email medical@medsurg.com.au www.medsurg.com.au

Phone 1300medsurg (1300 633 787)
Fax 07 3859 2995

