



Rapid, Simple and Reliable

Low-Sensitivity Urine Pregnancy Tests

Tests for the confirmation of Successful Medical Termination of Pregnancy

AccuPreg-Screen®

Self-Testing. This is not a classical pregnancy test.

Your healthcare professional has provided this brochure as it contains important information about the AccuPreg-Screen®.

This self-test is designed to help determine the outcome of a first-trimester Termination of Pregnancy (ToP) or miscarriage by measuring hCG levels above 1000mIU/mL.

The AccuPreg-Screen® is a rapid, simple, and reliable manual qualitative immunochromatographic screening test that, under the guidance of your clinician, can be performed in the comfort of your own home.

To obtain the Basin Medical AccuPreg-Screen® test, please consult your clinician.

The AccuPreg-Screen® test is:

-  Simple and easy to use.
-  Rapid, works in 5 to 10 minutes.
-  A low-sensitivity urine test.
-  Can be done at your discretion, either at home or wherever is most convenient for you.

How does the AccuPreg-Screen® test work?

When a woman becomes pregnant, her body starts producing a hormone called human Chorionic Gonadotropin (hCG), which is excreted in urine. The concentration of hCG in the urine increases over time. The AccuPreg-Screen® test utilises a pair of antibodies specifically designed to detect the hCG hormone in urine. This test checks whether the hCG level is 1000 mIU/mL or higher.

If a treatment for termination of pregnancy or a miscarriage is successful, the hCG level should decrease to below 1000 mIU/mL within 10 to 20 days. It's important to note that AccuPreg-Screen® will not detect hCG levels below 1000 mIU/mL, while traditional pregnancy tests may still show a positive result at those lower levels. Even if the test is negative, you should still contact your health care provider if you have symptoms of pregnancy, such as scanty bleeding, no period, tender breasts, feeling sick, and a growing tummy.




How to use the AccuPreg-Screen® test

The AccuPreg-Screen® test box contains one (1) hermetically sealed aluminium pouch containing:


- 1 x AccuPreg-Screen test
- 1 x desiccant bag
- 1 x Instructions for use leaflet



- 1 Remove the test from the pouch
- 2 Remove the cap from the test
- 3 Place the absorbent tip directly under your urine stream for 5 to 10 seconds.
- 4 Replace the cap on the test.
- 5 Wait 5 to 10 minutes.
- 6 Read the results. *Do not interpret results after 10 minutes.*

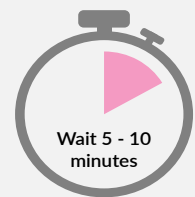

Results Indicator

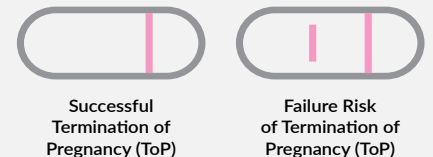

AccuPreg-Screen® is manufactured in France by Veda Lab
Rue De L'Expansion Zat Du Londeau Cerise, BP 181 Alencon, Cedex, 61006 France
AccuPreg-Screen® TGA Registration ARTG Entry: 490269



Medical and Surgical Requisites Pty Ltd

Call within Brisbane 07 3859 2900 33 Fulcrum Street Richlands Queensland 4077 www.medsurg.com.au
Phone 1300medsurg (1300 633 787) Email medical@medsurg.com.au Fax 07 3859 2995

