


# NEW! HEINE® iC 1

## 21<sup>st</sup> Century Digital Dermatoscopy

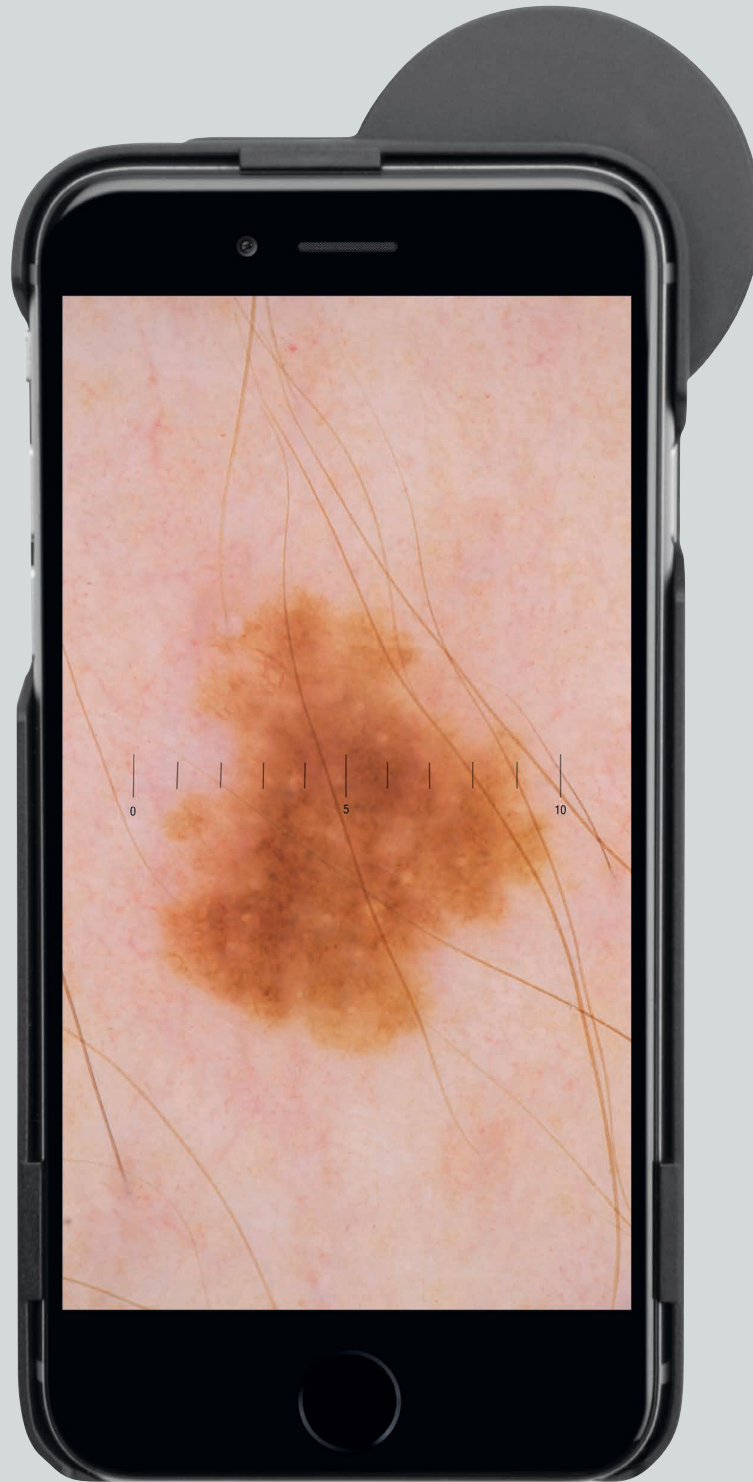

Established 1922

**Medical and Surgical Requisites Pty Ltd**

Call within Brisbane 07 3859 2900 33 Fulcrum Street Richlands Queensland 4077 [www.medsurg.com.au](http://www.medsurg.com.au)  
Phone 1300medsurg (1300 633 787) Email [medical@medsurg.com.au](mailto:medical@medsurg.com.au) Fax 07 3859 2995



# NEW! HEINE® iC 1 – 21<sup>st</sup> Century Digital Dermatoscopy

The HEINE iC 1 is a new image capturing device to record clinical and dermatoscopic images for documenting, monitoring and supporting differential diagnosis of melanoma and other skin lesions. Innovative HEINE LED<sup>HQ</sup> illumination combined with precision optics provide high resolution and reflex-free images with an Apple iPhone 6® or iPhone 6s®. The integrated Toggle Function allows you to easily alternate between polarised and non-polarised views.



Future compatibility – with exchangeable mobile phones cases

See, or show your patients examination results directly on the screen.

Onboard App  
Download on the App Store

Latest Li-Polymer battery technology with USB charging function and charge status indicator

Changeable and cleanable contact plate made of high quality glass. With scale for exact image comparisons over time.

Up to 40x magnification with HEINE precision optics and up to 12 Megapixel image quality (with iPhone 6s<sup>®</sup>\*)

Homogeneous LED<sup>HQ</sup> illumination

Toggle Function to conveniently alternate between polarised images and immersion images

[05]

## [ ACCESSORIES ]

<b>USB cord with E4-USB Medical approved plug-in power supply</b> [01] + [02] Medical approved plug-in power supply with 5V and 500mA; changeable plug adapter	Art.Nr: X-000.99.303
<b>Contact plate</b> with scale [03]	Art.Nr: K-000.34.201
<b>Mobile phone case</b> , currently available for Apple® iPhone 6® and iPhone 6s® [04]	Art.Nr: K-000.34.250



## [ TECHNICAL DATA ]

Catalogue number	<b>K-270.28.305</b> <b>iC 1 set/6</b> [05] complete with: HEINE iC 1 with contact plate with scale, mobile phone case iC 1/6 currently available for Apple® iPhone 6® and iPhone 6s®, USB cord with medical approved plug-in power supply
GTIN	4053755186044
Available versions	iC 1/6 for iPhone 6, iPhone 6s
Weight	140g
Dimensions	140x65x42 mm up to 155x75x50 mm
Users manual	DE, GB, FR, ES, IT, SE, NL, DK, FI, PT
Operating conditions	+10 °C to +35 °C, 30 % to 75 % RH, 700hPa to 1060hPa
Terms of transport	-20 °C to +50 °C, 45 % to 85 % RH, 500hPa to 1060hPa
Electrical data	micro USB (5V) <500 mA, charging time typ. 2,5h, operating time typ. 3,5h Li-Polymer rechargeable battery
Optical data	Optical magnification = 15x, digital magnification = 40x (with iPhone 6® and 24" display), 21 mm field of view
LED	4x high power LEDs, colour temp. typ. 5000K, working life typ. 100000h, typ. CRI ≥80

\*The iPhone 6® or iPhone 6s® is not included in the scope of delivery.  
iPhone is a trademark of Apple Inc., registered in the U.S. and other countries



## Medical and Surgical Requisites Pty Ltd

Call within Brisbane 07 3859 2900 33 Fulcrum Street Richlands Queensland 4077 www.medsurg.com.au  
Phone 1300medsurg (1300 633 787) Email medical@medsurg.com.au Fax 07 3859 2995

