


HEINE QUALITY
MADE IN GERMANY

HEINE® iC 2 Funduscope
Digital Ophthalmoscopy

NEW!

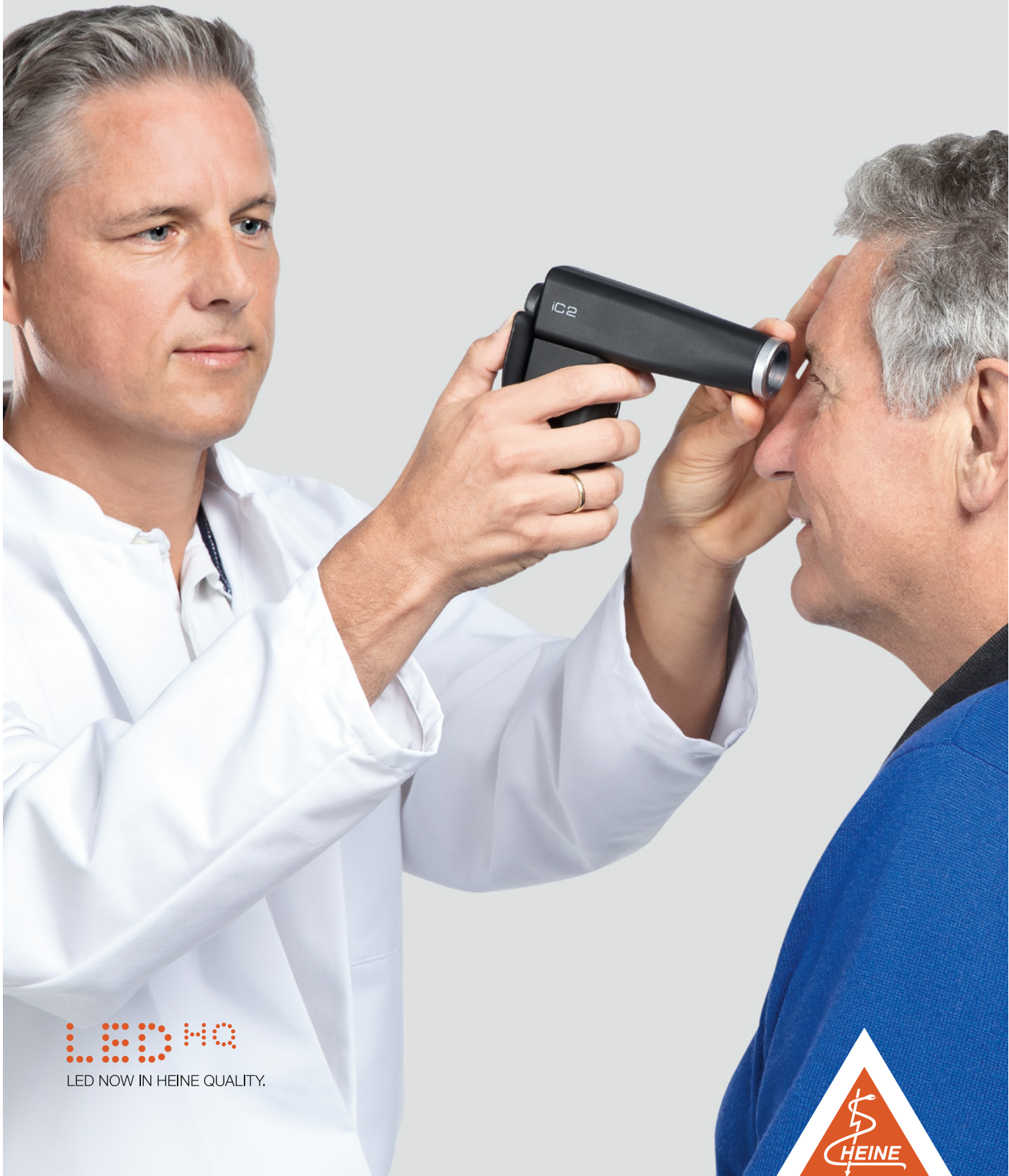

LED HO

LED NOW IN HEINE QUALITY.



Established 1922

Medical and Surgical Requisites Pty Ltd

Call within Brisbane 07 3859 2900
Phone 1300medsurg (1300 633 787)

33 Fulcrum Street Richlands Queensland 4077
Email medical@medsurg.com.au

www.medsurg.com.au
Fax 07 3859 2995



The new HEINE iC2 Funduscope is the next generation of direct, digital Ophthalmoscopes. Conveniently capture high quality fundus images.




- ✓ **Mobile and light**
Weighing only 300g, compact in size and supporting wireless operation, this instrument is ideal for convenient use even for patients who are confined to bed or have a disability.
- ✓ **High-quality fundus images**
Early identification of diseases of the eye. Ideal also for consultations. Documentation of current condition and development over time. Excellent image quality with up to 12 MP**.
- ✓ **Large field of view**
Field of view up to 34° (approx. 4x as large as the conventional direct ophthalmoscope with typically 9°).
- ✓ **Compensation of patient's defective vision**
Stepless manual dioptric adjustment of +/-15 dpt. and auto focus with +/-3 dpt.
- ✓ **Easy imaging of the periphery**
Scan the periphery as with a direct ophthalmoscope.
- ✓ **Photobiological Protection for the patient's eyes**
Complete brightness control and LED illumination which meets the highest demands (ISO 15004-2, Group I, international standard for the light hazard protection).
- ✓ **Sharp, shadow-free images**
Fine-tuned interaction with HEINE LED^{HQ} illumination with precision optics.




LED NOW IN HEINE QUALITY.


Examples of fundus images with the HEINE iC2 Funduscope

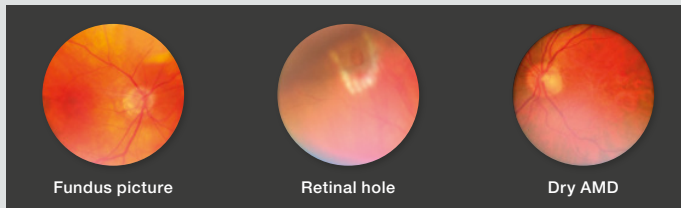

- ✓ **iC2 App**
Simple mapping of images. E-mail and zoom functions. Storage of additional patient data. Always up to date thanks to updates from the App Store*.

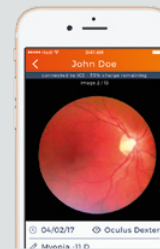

High-resolution images

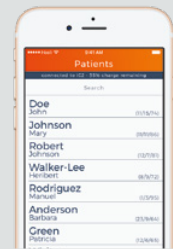

Select existing patients



Explain findings to patients

Ordering information:

| | |
|---|---------------------|
| iC2 Set/7 for iPhone 7 and 8 | |
| iC2 with mounting case smartphone iC 2/7 for Apple® iPhone® 7 and 8*, USB cord with E4-USB medical approved plug-in power supply and case | C-152.40.303 |
| iC2 Set/6 for iPhone 6/6s | |
| See above, but with mounting case smartphone iC 2/6 for Apple iPhone 6/6s* | C-151.40.303 |
| iC2 Set/5 for iPhone 5s and SE | |
| See above, but with mounting case smartphone iC 2/5 for Apple iPhone 5s and SE* | C-150.40.303 |

Due to local regulations, some products may not be available in your country. For additional product configurations and further HEINE Quality Products please visit www.heine.com

* The Apple iPhone is not included in the scope of delivery. Apple, the Apple logo, iPhone and iPod are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

** With iPhone from Version 6s

HEINE Optotechnik GmbH & Co. KG, Kientalstr. 7, 82211 Herrsching, Germany
Tel. +49 (0) 8152/38-0, Fax +49 (0) 8152/38-202, E-Mail: info@heine.com, www.heine.com



We reserve the right to change specification without notice.
© 2018 HEINE Optotechnik GmbH & Co. KG. All rights reserved. A-000.00.219 e 06/18.



Medical and Surgical Requisites Pty Ltd

Call within Brisbane 07 3859 2900 33 Fulcrum Street Richlands Queensland 4077 www.medsurg.com.au
Phone 1300medsurg (1300 633 787) Email medical@medsurg.com.au Fax 07 3859 2995

